

# SENTINEL™ feeding pump

The *Sentinel™* is a time-saving, easy-to-use enteral feeding pump that can be used in the administration of all enteral feeding formulas. The pump is easy to program, and many of its functions can be activated with a single touch. Reliable and accurate, the *Sentinel* is designed to give the user easy access to its many useful features.

◀ easy access ▶ ◀ easy-to-use ▶ ◀ easy-to-program ▶ ◀ reliable and accurate ▶

## Features and Benefits

- **Advanced Software** - pump odometer records and analyzes data for usage history
- **Battery Life** - internal monitor for battery life estimation
- **Large LED Display** - sizeable area for easy viewing
- **Alarm Identification** - differences in tone sequence make alarms easy to identify



## S<sup>2</sup> Feeding System™


*Curved over hook prevents dislodging tubing*

*Positive chamber retention holds Drip Chamber snug*

*Captive tube seating retains peristaltic section*

*Clover-shaped rotor for quick-see monitoring*

The Sentinel enteral feeding pump employs a patented system of mechanical features, known as the *S<sup>2</sup> Feeding System*, to ensure that the enteral pump spike set does not dislodge. This system works to protect against unintentional set removal.

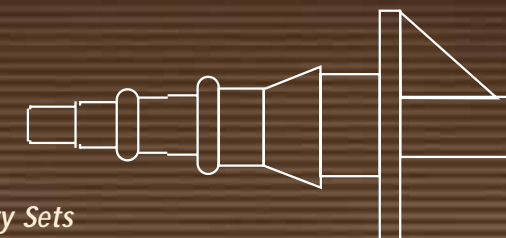

## Disposable Delivery Sets



## Bag Set



## Spike Set



**Secure Lock Tip**  
Interlocking tip design minimizes possibility of disconnection

Patent Pending

## 3 Step easy-to-use operation


1. Press **ON**.
2. Press **UP/DOWN** to set flow rate.
3. Press **RUN/STOP** to start/stop feeding.

**OST Medical**

**Impact resist-  
ant housing**  
withstands hard  
impacts

**Integral  
pole clamp**  
easy grip design

**Cap holder**  
for convenience  
and helps prevent  
contamination

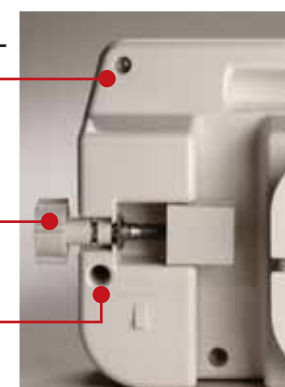

Rates	5-295 ml/hr in 1 ml increments
Dose Limit	1-2000 ml(s) in 1 ml increments/2000-9999 mls in 5 ml increments
Accuracy	+/- 10%
Battery	12 hours at 125 ml/hr
Recharge	Continuous while unit is plugged into AC outlet
Memory	Retained without power required
Size	7.5" x 9.1" x 5.2"
Weight	5.2 lbs.

### The Sentinel Enteral Feeding Pump from OST Medical

- Dose Check keeps a record of the accumulated amount of formula delivered over several feedings.
- Clear, Large Display Windows can be read from a distance. Displays are adjustable for brightness for day or night operation.
- The Clear Button clears all displays, except total accumulated dose, with one, four-second touch.
- Audible and Visible Alarms are readily differentiated. User is alerted when battery is low, when the formula container is empty, or the line is occluded, when the dose is complete, and in the event of formula excess flow.
- The Smooth, Plastic Case is easily cleaned with mild soap and water.
- Built-in Pole Clamp allows patient mobility.
- Pump is easily loaded with one hand.
- Memory is not erased in a power failure or when unit is turned off.
- Dose Limit can be set in 1 ml increments up to 2000 mls, then 5 ml increments up to 9999 mls.
- Programming Directions are located directly on the pump.
- Special Cap Holder for administration set on back of pump, keeps caps convenient and helps prevent contamination.

### Available Feeding Tubes

Part No.	Description	Packaging
E-1001	Spike Set	30 each
E-1200	1200 ml Bag Set with attached delivery set	30 each
E-2012	36" NG Tube, 12Fr., Non-Styleted, w/o Bolus	10 each
E-2014	36" NG Tube, 14Fr., Non-Styleted, w/o Bolus	10 each
E-1006	36" NG Tube, 6Fr., with Stylet and Bolus	10 each
E-1008	45" NG Tube, 8Fr., with Stylet and Bolus	10 each
E-1010	45" NG Tube, 10Fr., with Stylet and Bolus	10 each
E-1012	45" NG Tube, 12Fr., with Stylet and Bolus	10 each
E-1014	45" NG Tube, 14Fr., with Stylet and Bolus	10 each
E-1002	Y-Adapter	10 each

**SENTINEL™**  
feeding pump



*Simplicity, convenience  
and economy for*

*enteral feeding*

### OST Medical

OST Medical, Inc.  
11 Knight Street, Building F-23  
Warwick, RI 02886  
Tel.: 401-737-3774  
FAX: 401-737-4519

ISO 13485  
ISO 9001